


PALFIQUE ASTERIA

Excellent Physical Properties

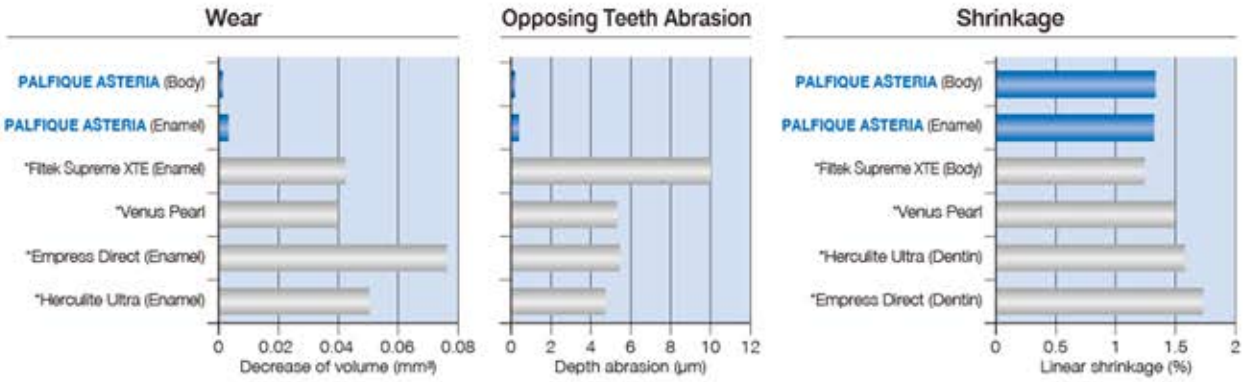

*Not a registered trademark of Tokuyama Dental
Source : Tokuyama Dental R&D

Product Packages Available



PALFIQUE ASTERIA SYRINGE ESSENTIAL KIT

Contents ..... 7 Syringes (2.1mL (4.0g) each)
Body; A1B, A2B, A3B, A3.5B, A4B
Enamel; NE, OcE
Clinical Guide



PALFIQUE ASTERIA SYRINGE Refill

Contents ..... 1 Syringe (2.1mL (4.0g))
Available shades: Body ..... A1B, A2B, A3B, A3.5B, A4B, B3B, BL
Enamel ..... NE, OcE, WE, YE, TE



PALFIQUE ASTERIA CUSTOM SHADE GUIDE

Contents ..... 15 Shade Guide Handles, Mold, Holder

Custom Shade Guide allows you to build and personalize according to your own preferences.
Customize your palette and create the exact shade of your composite.

ESTELITE COLOR KIT (Light-cured Characterizing Material)

Contents ..... 13 Syringes (0.5mL (0.9g) each)
15 Syringe Tips and Caps, 5 Disposable Dishes
1 Light Shield Cover.

ESTELITE COLOR Refill

Contents ..... 1 Syringe (0.5mL (0.9g))
2 Syringe tips and caps
Available shades ..... Clear, White, Blue, Yellow, Gray, Ochre, Dark Brown, Red, Lavender, Pink Opaque, High Chroma Opaque, Medium Chroma Opaque, Low Chroma Opaque

Brush No.24

Contents ..... 5 Double-Ended Brushes

Specifications are subject to change without obligation on the part of the Manufacturer.

Manufactured by
Tokuyama Dental Corporation

38-9, Taitou 1-chome, Taitou-ku, Tokyo, Japan
TEL: +81-3-3835-2261 FAX: +81-3-3835-2265
URL: www.tokuyama-dental.com

190523



SIMPLY AESTHETIC


PALFIQUE ASTERIA universal composite system is focused on simplified restorative therapy and outstanding aesthetic results.

Invisible Margins

PALFIQUE ASTERIA introduces a new simplified 2-step layering concept. This comprehensive system is comprised of minimal shades; just 7 Body shades and 5 Enamel shades cover the entire dental shade range.

BODY Shades

The Body shades have excellent blending ability with less width of margin bevel thanks to their state-of-the-art optical properties. The Body shades provide some translucency with sufficient opacity to avoid shining through without the use of opaque or dentin shades.
A1B-A4B shades blend with most natural dentition. The supplemental Body shade BL is recommended for high value bleached teeth and B3B is for yellowish teeth.

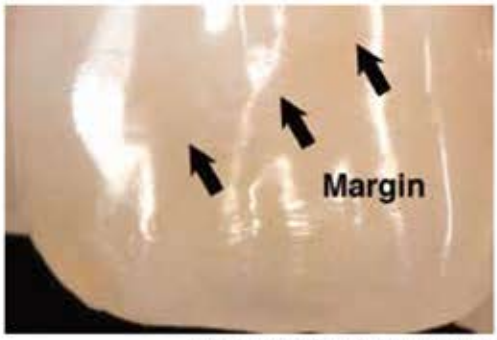

Case courtesy of Dr. Noboru Takahashi




A1B A2B A3B A3.5B A4B



BL B3B

ENAMEL Shades

The Enamel shades have compatible translucency as a substitute for natural enamel. The primary use of NE shade is for the anterior teeth and OcE is for the posterior occlusal area. The supplemental Body shade TE is for highly translucent anterior teeth, WE is for whitish enamel and YE is for discolored orangish enamel.

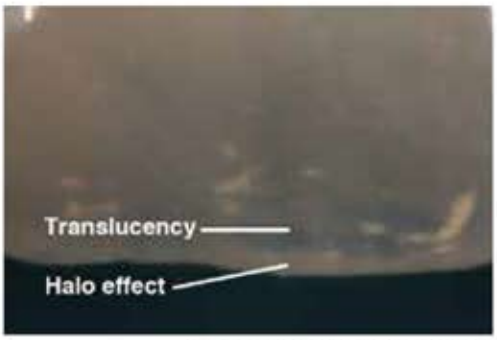

Case courtesy of Dr. Noboru Takahashi



NE OcE TE WE YE

Benefits

- **Outstanding Aesthetics**
  Invisible Margins
  Superior & Long Lasting Gloss Retention
- **Simplified Layering**
- **Optimal Handling**
  Non-Stick & Easy to Sculpt
  Ample Working Time
  (90sec. under ambient light (10,000 lux))






Simplified Layering

2-Step Layering Concept

Class I

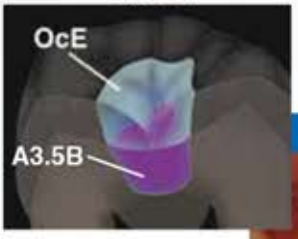


Postoperative view



Case courtesy of Dr. Noboru Takahashi

Class IV

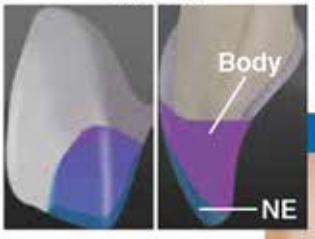

Postoperative view




Case courtesy of Dr. Noboru Takahashi

Superior & Long Lasting Gloss

The Supra-Nano Spherical filler in PALFIQUE ASTERIA facilitates a very smooth surface with superior gloss that is easily obtained with polishing and is sustained long term.


Filler SEM Images

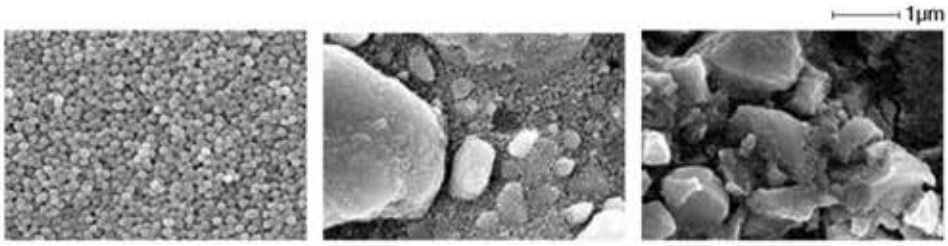

PALFIQUE ASTERIA | *Filtek Supreme XTE | *Venus Pearl

Gloss Retention

(Thermal cycle 10,000 times)

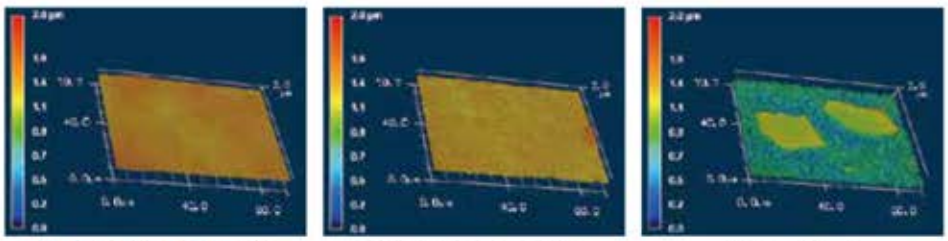

PALFIQUE ASTERIA | *Filtek Supreme XTE | *Venus Pearl

*Not a registered trademark of Tokuyama Dental

Source: Tokuyama Dental R&D